REVANESSE® SHAPE+™ with Lidocaine Product Specifications

1-844-226-8277
info@medsupplysolutions.com
medsupplysolutions.com


REVANESSE® SHAPE™+

Instructions for Use
Gebrauchsanweisung
Mode d'Emploi
Istruzioni per l'Uso
Instrucciones de Uso
Instruções de Uso
Οδηγίες χρήσης
Инструкции за употреба
Instrucțiuni de Utilizare
Instrukcje Stosowania
Návod k Použití
Bruksanvisning
Instruksjoner for Bruk
Brugervejledning

Instructies voor het Gebruik
Használati útmutató
Kasutusjuhend
Käyttöohjeet
Lietošanas instrukcija
Naudojimo instrukcija
Istruzioni għall-Użu
Návod na Použitie
Navodila za Uporabo
使用説明書
Инструкция по
применению
Kullanma Kılavuzu
تعليمات الاستعمال

Table of Contents:

Assembly of Needle to Syringe.....	3
Canada & Rest of the World.....	8
English (en).....	8
Français (fr).....	10
Español (es).....	12
Русский (ru).....	14
中文 (zh).....	16
EU.....	18
English (en).....	18
Français (fr).....	20
Deutsch (de).....	22
Italiano (it).....	24
Español (es).....	26
Português (pt).....	28
Ελληνικά (el).....	30
български (bg).....	32
Româna (ro).....	34
Polski (pl).....	36
Český (cs).....	38
Svenska (sv).....	40
Norsk (no).....	42
Dansk (da).....	44
Nederlands (nl).....	46
Magyar (hu).....	48
Eesti (et).....	50
Suomi (fi).....	52
Latviešu (lv).....	54
Lietuvių (lt).....	56
Malti (mt).....	58
Slovenský (sk).....	60
Slovenščina (sl).....	62
Русский (ru).....	64
Türkçe (tr).....	66
عربي (ar).....	68
Symbols.....	70
Contact.....	74

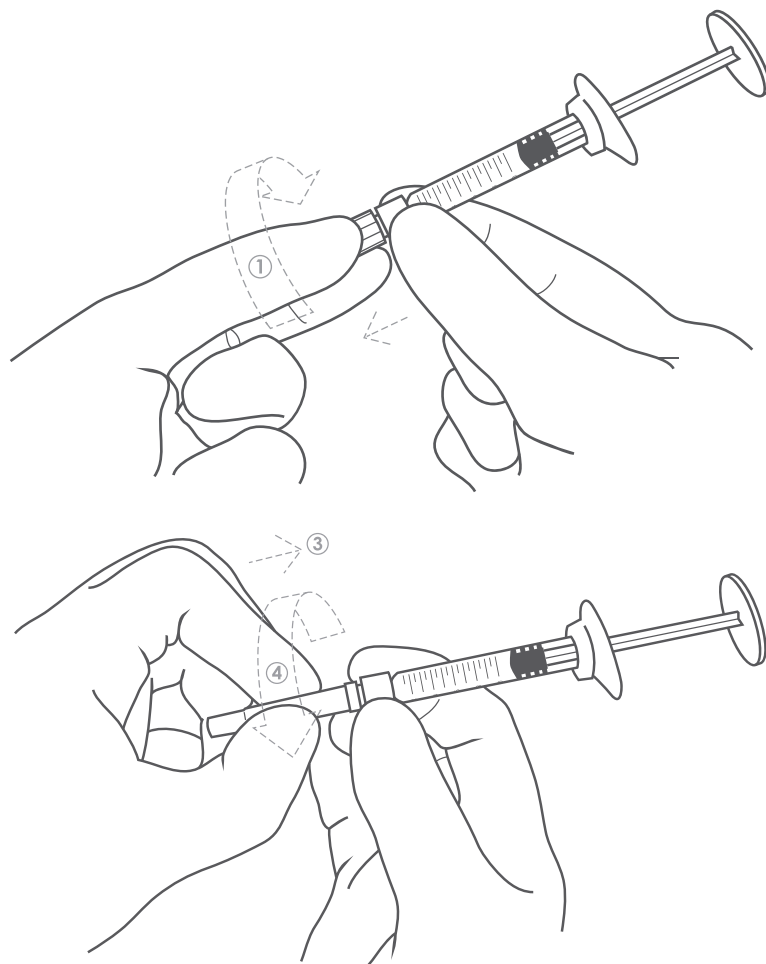ASSEMBLY OF NEEDLE TO SYRINGE:

1. Remove needle guard cap and unscrew the glass syringe end cap.
2. Firmly grip the glass syringe Luer-lock with thumb and forefinger.
3. Align the needle to the glass syringe Luer-lock.
4. Hold the glass syringe still and rotate the needle until it is tightly affixed.
5. Immediately prior to injection, remove needle guard by pulling forward (do not rotate).

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COMPOSITION

Cross-linked hyaluronic acid.....25mg/ml
 Lidocaine.....0.3% w/w
 In phosphate buffered saline
 [Cross linked with Butandiol-diglycidylether (BDDE)]

DESCRIPTION

Revanesse® Shape™+ is a colorless, odorless, transparent and aqueous gel of synthetic origin. The gel is stored in a pre-filled disposable syringe. Each box contains two 1.2ml syringes of Revanesse® Shape™+ along with up to two sterilized 27G needles.

APPLICATION RANGE / INDICATIONS

Revanesse® Shape™+ is a cross linked hyaluronic acid gel that is indicated for the treatment of medium to deep facial rhytides and lipoatrophy by injection into the skin. Implantation life is dependent on depth and location of injection, and averages 6-12 months.

ANTICIPATED SIDE EFFECTS

Physicians must inform patients that with every injection of Revanesse® Shape™+ there are potential adverse reactions that may be delayed or occur immediately after the injection. These include but are not limited to:

- Injection-related reactions might occur, such as transient erythema, swelling, pain, itching, discoloration or tenderness at the injection site. These reactions may last for one week.
- Nodules or induration are also possible at the injection site.
- Poor product performance due to improper injection technique.
- Glabellar necrosis, abscess formation, granulomas and hypersensitivity have all been reported with injections of hyaluronic acid products. It is important for physicians to take these reactions into account on a case by case basis.

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CONTRAINDICATIONS

- Contains lidocaine and is contraindicated for patients with a history of allergies to such material.
- Do not inject Revanesse® Shape™+ into eye contours (into the eye circle or eyelids).
- Pregnant women, or women during lactation should not be treated with Revanesse® Shape™+.
- Revanesse® Shape™+ is only intended for intradermal use and must not be injected into blood vessels. This may occlude and could cause an embolism.
- Patients who develop hypertrophic scarring should not be treated with Revanesse® Shape™+.
- Contains trace amounts of gram positive bacterial proteins, and is contraindicated for patients with a history of allergies to such material.
- Never use Revanesse® Shape™+ in conjunction with a laser, intense pulsed light, chemical peeling or dermabrasion treatments.
- People under the age of 18 should not be treated with Revanesse® Shape™+.
- Patients with acne and/or other inflammatory diseases of the skin should not be treated with Revanesse® Shape™+.
- Patients with unattainable expectations.
- Patients with auto-immune disorders or under immunotherapy.
- Patients with multiple severe allergies.
- Patients with acute or chronic skin disease in or near the injection sites.
- Coagulation defects or under anti-coagulation therapy.
- Patients with sensitivity to hyaluronic acid.

It is imperative that patients with adverse inflammatory reactions that persist for more than one week report this immediately to their physician. These conditions should be treated as appropriate (i.e: corticosteroids or antibiotics). All other types of adverse reactions should be reported directly to the authorized distributor of the Revanesse® family of products and / or to Prolenium Medical Technologies Inc. directly.

ADMINISTRATION & DOSAGE

- Revanesse® Shape™+ should only be injected by or under the direct supervision of qualified physicians who have been trained on the proper injection technique for filling facial wrinkles.
- Before patients are treated they should be informed of the indications of the device as well as its contraindications and potential undesirable side effects.
- The area to be treated must be thoroughly disinfected. Be sure to inject only under sterile conditions.
- Inject the product slowly and apply the least amount of pressure necessary.
- Revanesse® Shape™+ and needles packaged with it are for single use only. Do not re-use. If re-used, there is a risk of infection or transmission of blood born diseases.
- Keep the product at room temperature for 30 minutes prior to injection.
- If the skin turns a white color (blanching), the injection should be stopped immediately and the area should be massaged until the skin returns to its normal color.
- Before injecting, press on the plunger of the syringe until a small drop is visible at the tip of the needle.

PRECAUTIONS

- Revanesse® Shape™+ should not be injected into an area that already contains another filler product as there is no available clinical data on possible reactions.
- Revanesse® Shape™+ should not be injected into an area where there is a permanent filler or implant.
- Hyaluronic acid products have a known incompatibility with quaternary ammonium salts such as benzalkonium chloride. Please ensure that Revanesse® Shape™+ never comes into contact with this substance or medical instrumentation that has come into contact with this substance.
- Revanesse® Shape™+ should never be used for breast enlargement, or for implantation into bone, tendon, ligament or muscle.
- Avoid touching the treated area for 12 hours after injection and avoid prolonged exposure to sunlight, UV, as well as extreme cold and heat.
- Until the initial swelling and redness have resolved, do not expose the treated area to intense heat (e.g. solarium and sunbathing) or extreme cold.
- If you have previously suffered from facial cold sores, there is a risk that the needle punctures could contribute to another eruption of cold sores.
- If you are using aspirin, non-steroidal anti-inflammatory medications, St. John's Wort or High doses of Vitamin E supplements prior to treatment or any similar medications be aware that these may increase bruising and bleeding at the injection site.
- Based on a toxicological risk assessment for lidocaine, patients should be limited to 20 mL per 60 kg (130 lbs) body mass per annum. The safety of injecting greater amounts has not been established.
- The safety for use in patients under 18 years or over 65 years has not been established.
- Patients who are visibly ill, with bacterial or viral infections, influenza, or active fever should not be treated until a resolution of their symptoms.
- In order to minimize the risks of potential complications, this product should only be used by health care practitioners who have appropriate training, experience, and who are knowledgeable about the anatomy at and around the site of injection.
- Health care practitioners are encouraged to discuss all potential risks of soft tissue injection with their patients prior to treatment and ensure that patients are aware of signs and symptoms of potential complications.

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- Introduction of product into the vasculature may lead to embolization, occlusion of the vessels, ischemia, or infarction.
- Rare but serious adverse events associated with the intravascular injection of soft tissue fillers in the face have been reported and include temporary or permanent vision impairment, blindness, cerebral ischemia or cerebral hemorrhage, leading to stroke, skin necrosis, and damage to underlying facial structures.
- Immediately stop the injection if a patient exhibits any of the following symptoms, including changes in vision, signs of a stroke, blanching of the skin, or unusual pain during or shortly after the procedure.
- Patients should receive prompt medical attention and possibly evaluation by an appropriate health care practitioner specialist should an intravascular injection occur.

The Revanesse® family of products should not be used in areas that have high vascularity. Use in these areas such as the glabella and nose region has resulted in cases of vascular embolization and symptoms consistent with ocular vessel occlusion (i.e.: blindness).

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NOTE: The correct injection technique is crucial to treatment success & patient satisfaction. Revanesse® Shape™+ should only be injected by a practitioner qualified according to local laws and standards.

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COMPOSITION

Cross-linked hyaluronic acid.....25mg/ml
 Lidocaine.....3mg/ml (0.3% w/w)
 In phosphate buffered saline
 [Cross-linked with Butanediol-diglycidylether (BDDE)]

DESCRIPTION

Revasse[®] Shape[™]+ is a colorless, odorless, transparent and aqueous gel of synthetic origin. The gel is stored in a pre-filled disposable syringe. Each box contains two 1.2ml syringes of Revasse[®] Shape[™]+ along with up to two sterilized needles.

APPLICATION RANGE / INDICATIONS

Application: Medium to deep facial rhytides.

Medical Indications: The products are space-occupying tissue reconstructive materials composed of a hyaluronic acid gel that is indicated for restoration of volume lost from lipoatrophy/lipodystrophy, and/or correction of contour deficiencies and anatomic deformities of either pathologic origin or after trauma, in soft tissue.

Intended patients are those desiring the correction of contour deficiencies and deformities in soft tissue, such as HIV-associated lipoatrophy and lipodystrophy.

Cosmetic Indication: Revasse[®] Shape[™]+ is indicated for the treatment of facial rhytides, volume restoration, lip augmentation, skin hydration and contouring of depressions by injection into tissue.

ANTICIPATED SIDE EFFECTS

Physicians must inform patients that with every injection of Revasse[®] Shape[™]+ there are potential adverse reactions that may be delayed or occur immediately after the injection. These include but are not limited to:

- Injection-related reactions might occur, such as transient erythema, swelling, pain, itching, discoloration or tenderness at the injection site. These reactions may last for one week.
- Nodules or induration are also possible at the injection site.
- Poor product performance due to improper injection technique.
- Glabellar necrosis, abscess formation, granulomas and hypersensitivity have all been reported with injections of hyaluronic acid products. It is important for physicians to take these reactions into account on a case-by-case basis.

Reactions thought to be of hypersensitivity in nature have been reported in less than one in every 1500 treatments. These have consisted of prolonged erythema, swelling and induration at the implant site.

These reactions have started either shortly after injection or after a delay of 2 – 4 weeks and have been described as mild or moderate, with an average duration of 2 weeks. Typically, this reaction is self-limiting and resolves spontaneously with time. However, it is imperative that patients with hypersensitivity type reactions contact their physician immediately for assessment. Patients with multiple allergic reactions should be excluded from the treatment.

CONTRAINDICATIONS

- Products are contraindicated in cases of hypersensitivity to all local amide-type anesthetics (not only to lidocaine as indicated in the package leaflet), as well as to all degrees of atrioventricular block, sinus node dysfunction, severe decline in myocardial contractility, porphyria, supraventricular tachycardias.
- Contains lidocaine and is contraindicated for patients with a history of allergies to such material.
- Do not inject Revasse[®] Shape[™]+ into eye contours (into the eye circle or eyelids).
- Pregnant women, or women during lactation should not be treated with Revasse[®] Shape[™]+
- Revasse[®] Shape[™]+ is only intended for intradermal use and must not be injected into blood vessels. This may occlude and could cause an embolism.
- Patients who develop hypertrophic scarring should not be treated with Revasse[®] Shape[™]+
- Contains trace amounts of gram-positive bacterial proteins and is contraindicated for patients with a history of allergies to such material.
- Never use Revasse[®] Shape[™]+ in conjunction with a laser, intense pulsed light, chemical peeling or dermabrasion treatments.
- People under the age of 18 should not be treated with Revasse[®] Shape[™]+
- Patients with acne and/or other inflammatory diseases of the skin should not be treated with Revasse[®] Shape[™]+
- Patients with unattainable expectations.
- Patients with auto-immune disorders or under immunotherapy.
- Patients with multiple severe allergies.
- Patients with acute or chronic skin disease in or near the injection sites.
- Coagulation defects or under anti-coagulation therapy.
- Patients with sensitivity to hyaluronic acid.

It is imperative that patients with adverse inflammatory reactions that persist for more than one week report this immediately to their physician. These conditions should be treated as appropriate (ie: corticosteroids or antibiotics). All other types of adverse reactions should be reported directly to the authorized distributor of the Revasse[®] family of products and / or to Prolenium Medical Technologies Inc. directly.

ADMINISTRATION & DOSAGE

- Revasse[®] Shape[™]+ should only be injected by or under the direct supervision of qualified physicians who have been trained on the proper injection technique for filling facial wrinkles.
- Before patients are treated, they should be informed of the indications of the device as well as its contraindications and potential undesirable side effects.
- The area to be treated must be thoroughly disinfected. Be sure to inject only under sterile conditions.
- Inject the product slowly and apply the least amount of pressure necessary.
- Revasse[®] Shape[™]+ and needles packaged with it are for single use only. Do not re-use. If re-used, there is a risk of infection or transmission of blood born diseases.
- Keep the product at room temperature for 30 minutes prior to injection.
- If the skin turns a white color (blanching), the injection should be stopped immediately, and the area should be massaged until the skin returns to its normal color.
- Before injecting, press on the plunger of the syringe until a small drop is visible at the tip of the needle.

PRECAUTIONS

- If intravascular injections are made by mistake, the toxic effect will be noticed within 1-3 minutes (so the doctor and the patient should be alert in case of error).
- The effect of lidocaine may be reduced if infusion occurs in an inflammation or infection area.
- Special attention should be given to patients with partial or complete atrioventricular block because local anesthetics may suppress myocardial conduction in patients with advanced liver disease or severe renal impairment, in patients with epilepsy, in patients with respiratory failure, in the elderly age, in patients with poor general health status, in patients receiving class III antiarrhythmics (eg amiodarone), who must be under close medical supervision including electrocardiogram due to the possible addition of cardiac effects, and finally in patients with acute porphyria.
- Lidocaine should be used with caution in patients receiving other topical anesthetics or agents with structural similarities to local amide-type anesthetics, e.g. certain antiarrhythmics, such as mexiletine and tocainide, since systemic toxic effects can be cumulative.
- In addition to the direct anesthetic effect, local anesthetics may exert a very mild effect on the cognitive function and movement coordination and may temporarily affect somatomotility and alertness. Depending on the dose, local anesthetics may have a very minor effect on mental function and may temporarily disturb movement and movement coordination.
- Revasse[®] Shape[™]+ should not be injected into an area that already contains another filler product as there is no available clinical data on possible reactions.
- Revasse[®] Shape[™]+ should not be injected into an area where there is a permanent filler or implant.
- Hyaluronic acid products have a known incompatibility with quaternary ammonium salts such as benzalkonium chloride. Please ensure that Revasse[®] Shape[™]+ never comes into contact with this substance or medical instrumentation that has come into contact with this substance.
- Revasse[®] Shape[™]+ should never be used for breast enlargement, or for implantation into bone, tendon, ligament or muscle.
- Avoid touching the treated area for 12 hours after injection and avoid prolonged exposure to sunlight, UV, as well as extreme cold and heat.
- Until the initial swelling and redness have resolved, do not expose the treated area to intense heat (e.g. solarium and sunbathing) or extreme cold.
- If you have previously suffered from facial cold sores, there is a risk that the needle punctures could contribute to another eruption of cold sores.
- If you are using aspirin, non-steroidal anti-inflammatory medications, St. John's Wort or High doses of Vitamin E supplements prior to treatment or any similar medications be aware that these may increase bruising and bleeding at the injection site.
- Based on a toxicological risk assessment for lidocaine, patients should be limited to 20 mL per 60 kg (130 lbs) body mass per annum. The safety of injecting greater amounts has not been established.
- The safety for use in patients under 18 years or over 65 years has not been established.
- Patients who are visibly ill, with bacterial or viral infections, influenza, or active fever should not be treated until a resolution of their symptoms.
- In order to minimize the risks of potential complications, this product should only be used by health care practitioners who have appropriate training, experience, and who are knowledgeable about the anatomy at and around the site of injection.
- Health care practitioners are encouraged to discuss all potential risks of soft tissue injection with their patients prior to treatment and ensure that patients are aware of signs and symptoms of potential complications.

WARNINGS

Confirm that the seal on the box has not been broken and sterility has not been compromised. Confirm that the product has not expired. Product is for single use only; do not re-use. If re-used, there is a risk of infection or transmission of blood born diseases.

- Introduction of product into the vasculature may lead to embolization, occlusion of the vessels, ischemia, or infarction.
- Rare but serious adverse events associated with the intravascular injection of soft tissue fillers in the face have been reported and include temporary or permanent vision impairment, blindness, cerebral ischemia or cerebral hemorrhage, leading to stroke, skin necrosis, and damage to underlying facial structures.
- Immediately stop the injection if a patient exhibit any of the following symptoms, including changes in vision, signs of a stroke, blanching of the skin, or unusual pain during or shortly after the procedure.
- Patients should receive prompt medical attention and possibly evaluation by an appropriate health care practitioner specialist should an intravascular injection occur.

The Revasse[®] family of products should not be used in areas that have high vascularity. Use in these areas such as the glabella and nose region has resulted in cases of vascular embolization and symptoms consistent with ocular vessel occlusion (i.e.: blindness).

SHELF LIFE & STORAGE

Expiry is indicated on each individual package. Store between 2°-25° C, and protect from direct sun light and freezing.

NOTE: The correct injection technique is crucial to treatment success & patient satisfaction. Revasse[®] Shape[™]+ should only be injected by a practitioner qualified according to local laws and standards.

The graduation on the syringe is not precise and should be used as a guide only. The amount of material to be injected is best determined by visual and tactile assessment by the user.

MANUFACTURER

Prolenium Medical Technologies, Inc.
 138 Industrial Parkway North, Aurora, ON
 L4G 4C3, Canada

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This document is a photographic image of the product insert provided for your convenience by Med Supply Solutions. Med Supply Solutions does not claim ownership of the content contained within this product insert. The intellectual property rights, including but not limited to copyrights, trademarks, and any other applicable rights, pertaining to the content of this product insert, remain the sole property of the original manufacturer. Med Supply Solutions accepts no liability for any inaccuracies, omissions, or discrepancies that may be present in this document. Your use of this product insert is subject to our terms and conditions, which can be found on our website. By accessing and utilizing this document, you acknowledge and accept the terms and conditions outlined therein.

1-844-226-8277
info@medsupplysolutions.com
medsupplysolutions.com

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Emergo Europe B.V.
MedEnvoy Global B.V.